

One-Piece Pouching System

Collect Supplies



Gather pouch, disposal bag, tissue and/or cloth for cleaning, and water. Optional: scissors, pen, measuring guide, accessories.

Remove Used Pouching System



Carefully lift an edge of the adhesive and peel downward, gently pushing the skin away from the skin barrier. Properly dispose of used pouch. Do not flush down toilet.

Clean and Dry Skin



Clean the skin with water or with gentle soap as recommended by your healthcare professional. Rinse and dry thoroughly. Do not use oils, powder, lotion, or lotion-type soaps on the skin.


Check Condition of Skin



Check to make sure the skin around your stoma (peristomal) is intact without irritation, rash or redness.*

*Contact your WOC/NSWOC nurse if you notice irritation, rash or redness

Prepare the Skin Barrier and Pouch

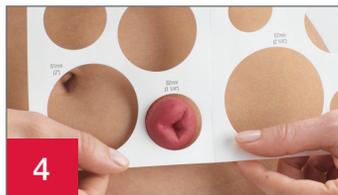

Measure the stoma using a stoma measuring guide before every barrier application (if using a pre-sized barrier opening, go to step 7).



Trace the determined barrier size opening onto the skin barrier. The skin barrier should fit where the skin and stoma meet. Verify that no skin is showing between the skin barrier and the stoma to help prevent leakage and skin irritation.



Using small curved scissors, carefully cut an opening in the skin barrier to match the pattern. Do not cut beyond the line on the release liner. Do not puncture the pouch.



Remove the release liner from the skin barrier.

Add Paste or Rings If Needed



If appropriate, apply a barrier ring or thin bead of paste around the opening following the instructions for use.

Apply the New Pouching System




Center the skin barrier opening over the stoma.



Apply gentle pressure to the skin barrier for about a minute for best adhesion.



For tape-bordered barriers, remove the backing paper on both sides and gently press adhesive against skin.

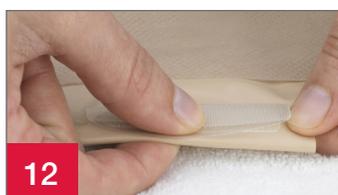

Close bottom of pouch. If using a drainable pouch, close the **Lock 'n Roll** microseal closure at the bottom of the pouch or apply a pouch clamp. Refer to the *Lock 'n Roll Microseal Closure Care Tip* for more information on how to use this closure system.



If using a urostomy pouch, make sure the drain valve is closed. The teardrop on the drain valve should face the body side of the pouch. If desired, apply the cap provided in the box.

One-Piece Pouching System

TIPS	LIST OF PRODUCTS	STOCK NO
<ul style="list-style-type: none">• The skin around your stoma should not show any signs of rash or sores• Irritated skin should be evaluated by a healthcare professional• If using a barrier ring or paste, apply it to the sticky side of the skin barrier next to the stoma opening as directed• Your stoma will probably be swollen after surgery. It may take several weeks or months for the stoma to shrink to its permanent size• The stoma may bleed slightly when you clean it• The use of skin protective wipes are not recommended under extended wear skin barriers as they may decrease wear time	One-Piece Pouching System	_____
	Accessories	
	Barrier rings	_____
	Ostomy belt	_____
	Paste	_____
	Lubricating deodorant	_____
Stoma powder	_____	
Other _____	_____	

Additional Instructions

For product questions, sampling needs, or detailed clinical questions concerning our products in the US, call **1.888.808.7456**. In Canada call **1.800.263.7400**.

Skin should be clean and dry before applying the skin barrier.

Prior to use, be sure to read the Instructions for Use for information regarding Intended Use, Contraindications, Warnings, Precautions, and Instructions.

The Hollister logo, Lock 'n Roll and "Healthy skin. Positive outcomes." are trademarks of Hollister Incorporated.
© 2018 Hollister Incorporated. 923069-1118 US-00667

Hollister Incorporated
2000 Hollister Drive
Libertyville, Illinois 60048
1.800.323.4060

www.hollister.com

Hollister Limited
95 Mary Street
Aurora, Ontario L4G 1G3
1.800.263.7400

